


Whittemore Veterinary Digital Platform Scale

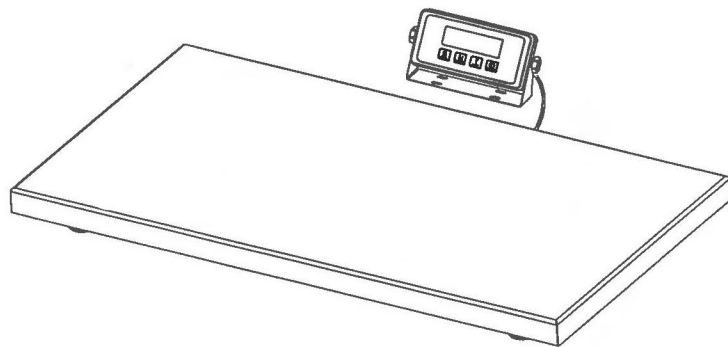

KEY DESCRIPTIONS:

HOLD	MODE	TARE	ON/OFF
------	------	------	--------

- 1. ON/OFF:**
Press the key to turn the scale on or off
- 2. MODE:**
Select weighing unit: kg and lb .
- 3. TARE:**
To deduct the weight of the container such as pan on the scale and the tared range is up to 100% of full capacity .
- 4. HOLD:**
Press this key to store weight value
- 5. POWER SWITCH :**
Press the switch to turn the power on or off which is on the back side of the indicator.

HOLD FUNCTION

- 1) Place the item on the Scale platform.
- 2) Wait for the stable weight to be displayed.
- 3) Press the HOLD key. The " √ " hold indicator will turn on.
- 4) Remove the item from the Scale platform. He item's weight reading will remain on the display. Now you can also select Unit (Pounds or Kilograms) by pressing [MODE] key.
- 5) Quit HOLD function: Once an item is weighed and the hold function is enabled, you want to quit HOLD function, press the HOLD key again to cancel the " √ " hold indicator and return the Scale to weight mode.


Calibration:

Press and hold [MODE] key then press [ON/OFF] key, the display will show CAL, then AD value, press [MODE] key when the stable indicator, the display will show flash calibration weight, place the shown weight on the platform, press [MODE] key when the stable indicator lights, the display will show AD value. Turn off the scale now, and calibration is finished.

FUNCTION SETTINGS

1. Press and hold [HOLD] key then press [ON/OFF] key to turn on the scale, the display will show A. oFF, press [MODE] then press [TARE] key to select auto power off time. press [HOLD] key to confirm
2. Press [HOLD] key again to enter next menu, the display will show b.L . press [MODE] then press [TARE] key to select auto backlight mode("On", "Off", "Auto" or "stb". press [HOLD] key to confirm
3. Press [HOLD] key again to enter next menu, The display will show Ani. Su, press [MODE] then press [TARE key to select ON (means animal weighing mode) or OFF (means normal weighing mode), press HOLD key to confirm.
4. Press [HOLD] key again to enter next menu, display will show FiLd.L press [MODE] then press [TARE] key to select select Filter rating among str, week. nor press HOLD key to confirm.
5. Press [HOLD] key again to enter next menu, the display will show LoCK.r press [MODE] then press [TARE] key to select animal scale anti-shock range among 10,30,50,oFF. press HOLD key to confirm.

Features:

1. Auto back light (selectable)
2. Auto shut off (selectable)
3. Hold function
4. Animal weighing
13. Display options: 6 digital LCD
14. Power: AC adapter 6V 100mA
15. Operating Temperature: 5~35°C
16. Operating Humidity: 25%~95%RH
17. Error Message indication:
 - Err-0 Overload
 - Err-5 Unstable when power on
 - Err-c Calibration code is not enough or operation error
 -  Low battery indication

Whittemore Enterprises, Inc
Sales, Service & Repair
(800) 999-2452
11149 Arrow Route
Rancho Cucamonga, CA 91730
www.wemed1.com
sales@wemed1.com
biomed@wemed1.com