


# Slender Arthroscope, size: 2.7 mm/3.0 mm

with **HOPKINS** Telescope. Diameter: 2.4 mm.  
For small joints.

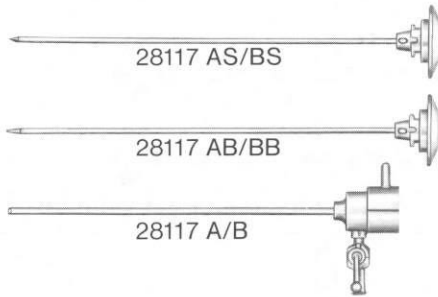


- 28116 A **HOPKINS** Straight Forward Telescope 0°. Fiber optic light transmission incorporated. Color code: green.
- 28117 A **Arthroscope Sheath, 2.7 mm Ø**, with one stopcock, with sharp obturator 28117 AS and blunt obturator 28117 AB, including one LUER-lock adaptor with stopcock 27502. Color code: green.
- 28116 B **HOPKINS** Forward-Oblique Telescope 30°. Fiber optic light transmission incorporated. Color code: red.
- 28117 B **Arthroscope Sheath, 2.7 mm Ø**, with one stopcock, with sharp obturator 28117 BS and blunt obturator 28117 BB, including one LUER-lock adaptor with stopcock 27502. Color code: red.

The 30° Forward-Oblique Telescope is the preferred telescope. It offers a straight forward as well as an oblique view. With rotation of this telescope you are afforded visualization in a direction opposite from that being observed, thereby increasing the range of vision.

The 30° slender arthroscope 28116 B/28117 B is very useful for the inspection of the posterior portion of the shoulder and elbow joints.

- Short-tipped obturators for safe puncture of thin walls of small joints.
- Ports to connect irrigation tubing and connection for the light cable are located in the same plane at a 6 o'clock position for the comfort of the surgeon.

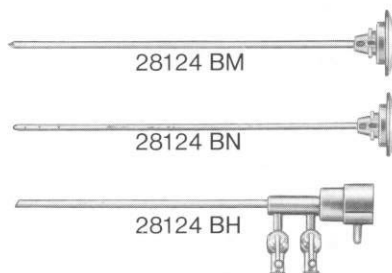

- 28124 BH **Arthroscope Sheath, 3 mm Ø**, with two stopcocks, with sharp obturator 28124 BM and blunt obturator 28124 BN, for use with 28116 B **HOPKINS** telescope 30°, including one LUER-lock adaptor with stopcock 27502 and one LUER-lock adaptor without stopcock 27500. Color code: red.
- 28117 Z Biopsy Forceps.

