

Simultaneous Lithotomy and Abduction!

Note: Stirrups are Sold in Pairs.

The Yellofin Stirrup represents the gold standard in lithotomy positioning. Yellofin Stirrups allow for easy adjustment of abduction and lithotomy while maintaining the sterile field. The boot design reduces pressure under the popliteal fossa and the superficial peroneal nerve. The Yellofin enables safe and easy positioning while providing enhanced surgical site access.

- **Boot Design** – Fin design helps reduce pressure under the popliteal fossa and superficial peroneal nerve.
- **Floating Boot** – Self-adjusting boot minimizes pressure on the calf, when moving the stirrup.
- **Lift-Assist** – Allows easy movement of the leg when placing it in the desired position.
- **Yellofin Reusable Pad** – Completely encapsulates the foot, ankle and calf for comfortable and secure positioning.
- **Squeeze Grip Handle** – Provides easy intra-operative adjustment without compromising the sterile field. Simply release the handle to secure the leg holder in all directions.
- **Lithotomy range** – Set the stirrup in any position between +90° to -35° lithotomy.
- **Abduction Range** – Set the stirrup in any position between +25° to -9° abduction.
- **Lithotomy & Length Indicators** – Visual indicators that help ensure precise positioning.

#O-YFES Yellofin Elite with Lift-Assist (500 lbs./227 kg Capacity)

#O-YFASI Yellofin with Lift-Assist (350 lbs./159 kg Capacity)

#O-YFSI Yellofin without Lift-Assist (350 lbs./159 kg Capacity)


Choose A Product	 Yellofin Elite with Lift-Assist #O-YFES 1 Pair	OR	 Yellofin with Lift-Assist #O-YFASI 1 Pair	OR	 Yellofin without Lift-Assist #O-YFSI 1 Pair
	 Yellofin Reusable Pad #O-AMA68S0.6S 1 Pair Required		 Easy Lock Blade Clamp #A-40040 (US), #A-40041 (Eur), #A-40042 (UK), #A-40043 (Jpn) 2 Required, Specify Rail Type		 Safety Drape #A-22000 Case of 10 Pairs
	 Stirrup Cart #A-30015		 Uro Catcher System #O-UC4 (US), #O-UC4E (Eur), #O-UC4UK (UK) (See Page 37 for More Info)		 Extended Warranty #A-XWARYFPRO Available in US Only



Easily Grip and Squeeze...



...Raise and Abduct...

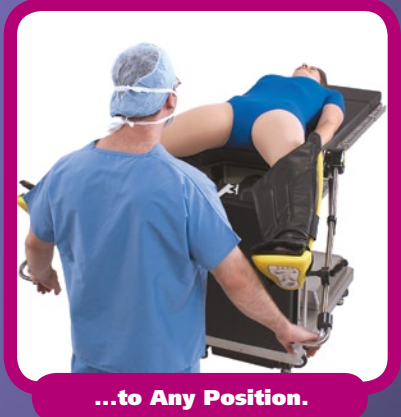

...to Any Position.



Squeeze and Position