


stryker[®]
ENDOSCOPY

THE SE5 SHAVER SYSTEM

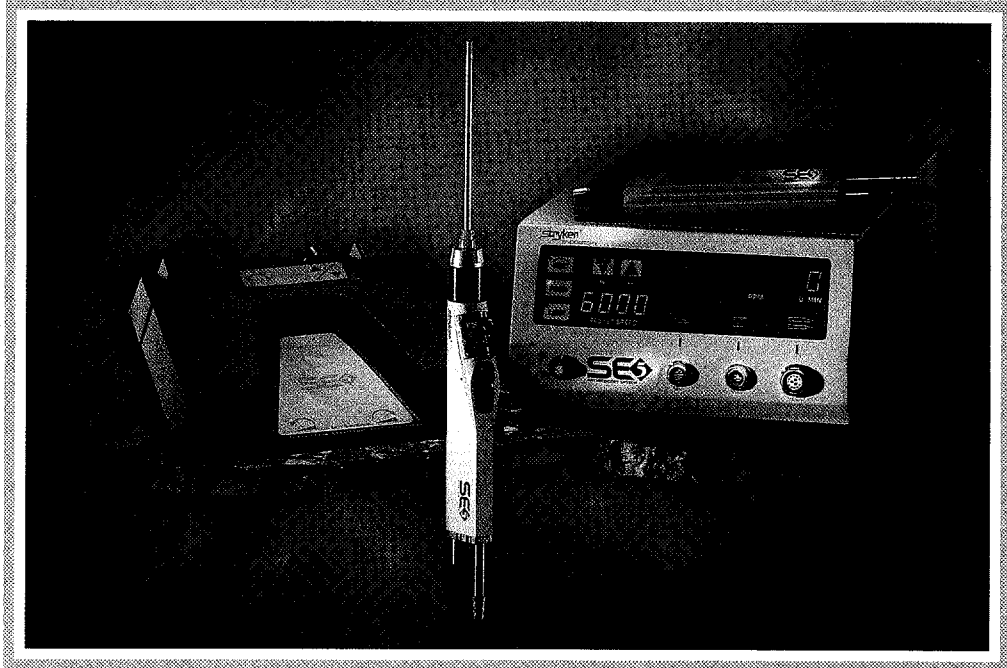

**THE SE5
SHAVER SYSTEM**


POWER
With 40 in.-oz. of torque and 6000 rpm, the SE5 Shaver offers you a tremendous combination of speed and power. The high speed and torque provide aggressive soft tissue resection and rapid burring. The shaver oscillates from 0 - 3000 rpm to provide quick soft tissue removal.

VERSATILITY
The SE5 Shaver offers you the ability to run 4 different handpieces and the ability to completely control your speed in the sterile field. The SE5 Console powers the SE5/QuadraCut Shaver Handpiece, the SE5 Hand-Controlled Shaver Handpiece, the SE5/QuadraCut Small Joint Handpiece, and the Bone Plug Handpiece.



The SE5/QuadraCut Shaver Handpiece is our traditional foot controlled handpiece and can be used in both the SE5 and QuadraCut ACL/Shaver System.

The SE5 Hand-Controlled Handpiece gives you the ability to control your handpiece with fingertip controls on the handpiece or by use of the footpedal.

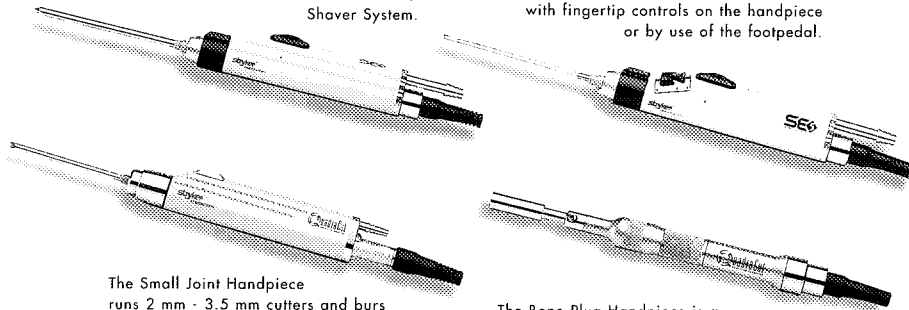

The Small Joint Handpiece runs 2 mm - 3.5 mm cutters and burs for your small joint arthroscopy needs.

The Bone Plug Handpiece is a circular oscillating saw for harvesting round bone plugs during patellar tendon harvesting when performing ACL reconstruction.

CONTROL

The Stryker Shaver Handpiece gives you the ability to completely control your speed easily and accurately with our variable footpedal. Simply press harder to go faster!



RELIABILITY

Surgeons have been relying on Stryker for quality products for over 40 years. The SE5 Shaver is our 5th generation shaver and continues our tradition of providing highly reliable products.

SE5 SPECIFICATIONS

SE5 Console	272-702
Size	9.5 in. (241.3 mm) width 14.1 in. (358.1 mm) depth 5.0 in. (127.0 mm) height
Weight	21.5 lbs.
Power Requirements	100-120V- 50-60 Hz 2.5A
SE5 Hand-Controlled Handpiece	272-704
Size	9.5 in. (241.3 mm) length 1.0 in. (25.4 mm) diameter
Weight	2.2 lbs. (1 Kg)
Speed	6,000 rpm
Torque	40 in.-oz.
QuadraCut Shaver Handpiece	275-701
Size	9.5 in. (241.3 mm) length 1.0 in. (25.4 mm) diameter
Weight	2.2 lbs. (1 Kg)
Speed	6,000 rpm
Torque	40 in.-oz.
QuadraCut Bone Plug Handpiece	275-705
Size	8.4 in. (213.4 mm) length 1.0 in. (25.4 mm) diameter
Weight	1.8 lbs. (0.82 Kg)
Speed	22,400 cpm
QuadraCut Small Joint Handpiece	275-601
Size	7.3 in. (185.4 mm) length 0.88 in. (22.4 mm) diameter
Weight	1.8 lbs. (0.82 Kg)
Speed	3,000 rpm
Torque	34 in.-oz
SE5 Footswitch	272-703
Size	8.87 in. (225 mm) width 7.25 in. (184 mm) depth 2.25 in. (57 mm) height
Weight	2.95 lbs. (1.34 Kg)
Operation Modes	forward reverse oscillate high/low