

OLYMPUS®

OLYMPUS ENF TYPE P3 RHINO-LARYNGOFIBERSCOPE

*Brightness And Resolution Enhancements
Assure Clear Observation For The Detection Of Disease
In Its Earliest Stages.*



Designed Specifically For Superior Observation And Diagnosis.



Olympus' improvement of the two crucial factors integral to fiberoptic utility assures greater accuracy and confidence for diagnosis or observation. Enhanced brightness and increased

resolution make the ENF TYPE P3 ideal for applications involving the nasal sinus, paranasal sinus, nasopharynx, larynx, and external ear.

Brightness increased approx. 48%*

A newly developed Olympus optical system provides an image field half again as bright as was previously possible. In addition to superior observation of the nasopharynx, the 85° field of view provides excellent visibility of the vocal cords.

Resolution improved approx. 33%*


By increasing the number of image guide fibers, Olympus has upgraded resolution by a third. This in turn allows the user to discern far smaller abnormalities, providing a much-needed boost in accuracy and precision for early diagnosis.

Angulation increased to 130° down

Olympus has added a full 30° of angulation capability to aid insertion and observation. Instead of

a conventional 100°, the ENF TYPE P3 can bend 130° down, making it especially useful for observing the nasopharynx and hypopharynx via nasal insertion.

*As compared to current model.

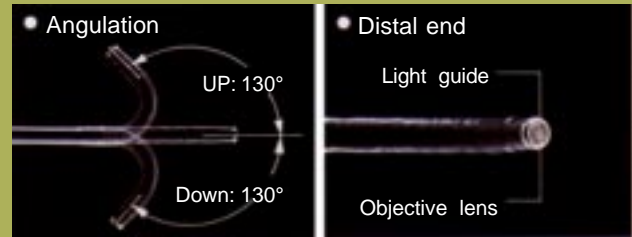

Specifications

Optical System	Field of view Direction of view Depth of field Illumination method	85° 0° (Forward viewing) 5~50mm Light guide system
Distal End	Outer diameter	3.4 mm
Bending Section	Range of tip bending	Up 130°, Down 130°
Insertion Tube	Outer diameter	3.6 mm
Working Length		255 mm
Total Length		485 mm
Photo Documentation	Still	35 mm SLR Camera SC35/OM-1 n (with OM Adapter SM-R)
	TV	OES TV System OLYMPUS OTV-S4

Specifications are subject to change without any obligation on the part of the manufacturer.



Photographic, Medical, Microscopic,
Industrial & Business Equipment

OLYMPUS

OLYMPUS OPTICAL CO., LTD.
San-Ei Building, 22-2, Nishi Shinjuku 1-chome, Shinjuku-ku, Tokyo, Japan
OLYMPUS OPTICAL CO.(EUROPA) GMBH.
Postfach 10 49 08 20034, Wendenstrasse 14-16 20097 Hamburg Germany
OLYMPUS AMERICA INC.
4 Nevada Drive, Lake Success, N Y 11042-1179. U S A
OLYMPUS KEYMED
KeyMed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England
OLYMPUS SINGAPORE PTE LTD
BLK 211, Henderson Road = 13-03, Henderson Industrial Park, Singapore 0315
OLYMPUS BEIJING REPRESENTATIVE OFFICE
Room No 3406 Beijing Jing Guang Centre, Hu Jia Lou, Chao Yang Qu, Beijing, China