

For the small-bone specialist.

POWERTOOLS

From the small-bone specialist.

SmartDriver

The versatility you want, the power you require for a wide range of procedures.

*Part of the
MicroAire SmartDrive™
System of power
instruments.*



6640 SMARTDRIVER

The 6640 SmartDriver with wire collet provides over 135 in. oz. of torque and over 1600 RPM of speed to insert K-wires up to .062" for fixation. Attachments are available for drilling, pinning, sawing and lavaging. Adapts automatically to wires and pins from .028" to .125". Variable speed trigger with "forward" and "reverse" switch on the handpiece permits one-hand operation. Choose electric or pneumatic.

MicroAire®
SURGICAL INSTRUMENTS

One powerful little company.

SmartDriver Attachments

HANDPIECE 6640

Speed: 1600 rpm nominal
Torque: Over 135 in. oz.
Cannulation: 3.2 mm (.125" or 1/8")

AUTOMATIC WIRE COLLET 6645

Automatic locking
Wire Advancement Tension Lever
Variable Speed Trigger
Accepts wires up to .062" (1.6mm)

AUTOMATIC PIN COLLET 6650

Automatic locking
Pin Advancement Tension Lever
Variable Speed Trigger
Accepts pins from .078" to .125"
(3.2mm or 1/8")

PULSE LAVAGE COUPLER 6640-474

Flow Rate: 500-850ml
Pulses: 5500 per minute
Speed: variable
Sterile, disposable tubing/nozzle sets available (See Pulse Lavage Brochure)

TRINKLE CHUCK 6660

Accepts drill bits, screwdrivers and other accessories with Trinkle fitting.

SAGITTAL SAW 6670

Speed: 18000 cpm
7° Blade Stroke
Saw Positions: 4
Blade Positions: 5

5/32" JACOBS CHUCK 6680


Accepts wires, pins, and drill bits up to 3.9mm (.156")

1/4" HIGH TORQUE JACOBS CHUCK 6685

Accepts wires, pins, and drill bits up to 6.4mm (.250")
Torque: 31 in. lb.
Speed: 400 rpm

SYNTHE CHUCK 6690

Accepts ASIF/AO (Synthes) twist drills.

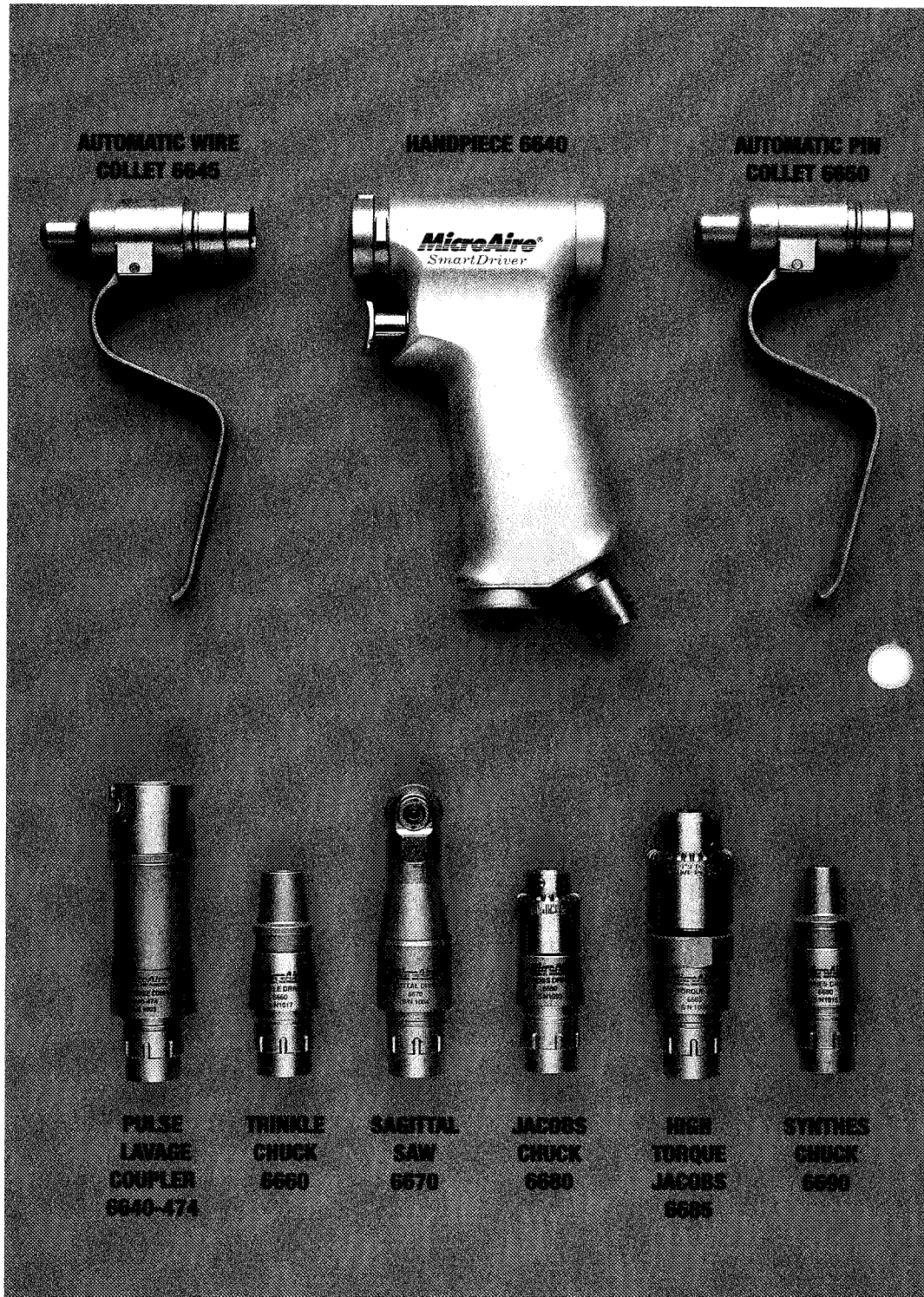

ORDERING
INFORMATION

CATALOG#	DESCRIPTION
6640	Pneumatic Handpiece (Choice of air hose connectors)
6640 E	Electric Handpiece
6640-474	Pulse Lavage Coupler (See Pulse Lavage Brochure for Accessories)
6645	Wire Collet (Accepts wires up to .062" or 1.6mm)
6650	Pin Collet (Accepts pins from .078" to .125" or 3.2mm)
6660	Trinkle Chuck

6670	Sagittal Saw
6680	5/32" Jacobs Chuck with key
6685	1/4" High Torque Jacobs Chuck with key
6690	Synthes® Chuck
6025	SmartDrive™ Electric Console
6006-000	Electric Handpiece Cord
9000-000	Air Hose for 6640-100