


USA Series™ Micro-H™ Operative Continuous Flow Hysteroscopes

A Simple Solution to a Difficult Problem—Micro Hysteroscopy

The CIRCON ACMI™ Micro-H Hysteroscopes offer a variety of configurations for diagnostic and therapeutic capabilities. The integrated small diameter hysteroscope (4.0 mm or 4.6 mm) utilizes advanced fiberoptics with diopter correction. Continuous flow is created by snapping on an outer sheath (4.7 mm or 5.5 mm). Large working channels (1.7 mm or 2.3 mm) and long working lengths (both 25 cm) make hysteroscopy simple and versatile in the office, surgicenter, or operating room. The working channels provide direct access for a variety of operative flexible instruments.

Micro-H™

GYH-4607 Hysteroscope with GYS-4607 Sheath

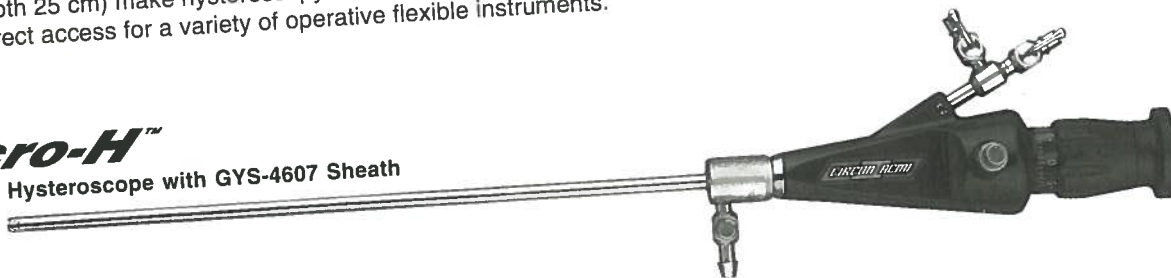

Features

- Small diameter: 4.0 mm or 4.6 mm
- Large working channel: 1.7 mm or 2.3 mm (5 Fr or 7 Fr)
- Continuous flow
- Advanced fiberoptic lens design
- Semi-rigid technology

Benefits

- Easier cervical insertion
- Provides access for flexible instruments
- Maintains clear field of view
- Produces sharp, clear images
- Proven durability

Hysteroscopes

GYH-4005 Micro-H Operative Continuous Flow Hysteroscope, 4.0 mm, 1.7 mm working channel

GYH-4607 Micro-H Operative Continuous Flow Hysteroscope, 4.6 mm, 2.3 mm working channel

Continuous Flow Outer Sheaths

The continuous flow outer sheaths snap over the Micro-H Hysteroscope and provide a means to drain distention medium.

GYS-4005 Micro-H Continuous Flow Outer Sheath, 4.7 mm, for GYH-4005

GYS-4607 Micro-H Continuous Flow Outer Sheath, 5.5 mm, for GYH-4607

SPECIFICATIONS

Working Length:	25 cm
Working Channel—	
GYH-4005:	1.7 mm / 5 Fr (GYH-4005)
GYH-4607:	2.3 mm / 7 Fr (GYH-4607)
Continuous Flow Sheath—	
GYH-4005:	4.7 mm (GYS-4005)
GYH-4607:	5.5 mm (GYS-4607)
Focusing Ocular:	Yes
Direction of View:	12°

COMPANION PRODUCTS

CS-B612	Biopsy Port Seals, self-sealing, 6–12 Fr (6/pkg)
GYN-6431	Flexible Gynecologic Biopsy Forceps, 6 Fr (2 mm), 43 cm
GYN-6432	Flexible Gynecologic Alligator Forceps, 6 Fr (2 mm), 43 cm
GYN-6433	Flexible Gynecologic Scissors, 6 Fr (2 mm), 43 cm
GYN-6436	Flexible Gynecologic Biopsy Forceps, 5 Fr (1.7 mm), 43 cm
GYN-6437	Flexible Gynecologic Alligator Forceps, 5 Fr (1.7 mm), 43 cm
GYN-6431S	Flexible Gynecologic Biopsy Forceps Stem Only, 6 Fr (2 mm), 43 cm
GYN-6432S	Flexible Gynecologic Alligator Forceps Stem Only, 6 Fr (2 mm), 43 cm
GYN-6433S	Flexible Gynecologic Scissors Stem Only, 6 Fr (2 mm), 43 cm
330H	Handle for GYN-6431S, -32S, -33S

