


Ritter 104

Examination table

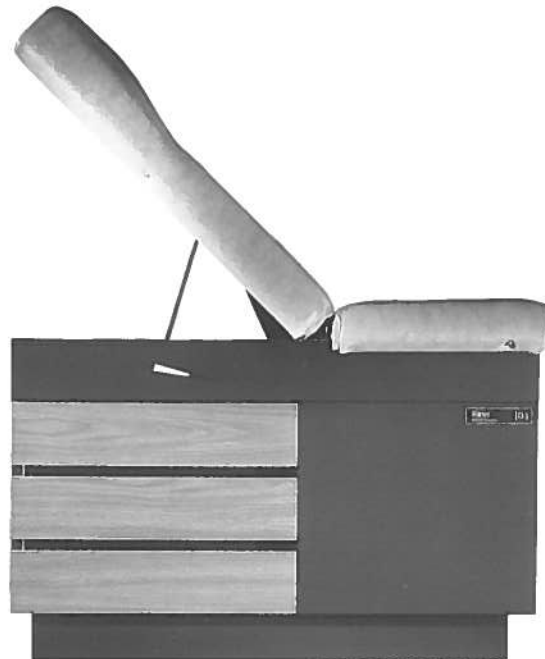

This is the most copied examination table in the world. And the competition still doesn't have it right. There are five big drawers to maximize available storage space, plus a pull-out step covered with Scotch-Tred® for safe mounting and dismounting.

Room at the top

The patients will like the comfort of the Soft-touch Top. The sculptured vinyl, stain-resistant top extends to 72" with a super-sturdy pull-out foot extension. Under the foot extension is a convenient pull-out treatment pan. The easily cleaned top has a full 1½" of foam padding with extra soft edges and corners.

Getting your back up

The manual counter-balanced back section positions easily



because the back is spring assisted. It can be adjusted quickly and easily to any position you choose. Under the back is a paper roll holder.

Attention to detail

The drawers use stainless steel roller bearing glides. The stops are positive and there's virtually no side-to-side play. The stirrups are easy to adjust, lock at any length in three lateral positions and store out of sight.

Options & Accessories

The combination Pelvic Tilt/Drawer Heater is a valuable factory-installed option. The heater is a blessing for physicians as well as patients while the pelvic tilt, set at 1°, can be positioned to 4° and 8° from either side of the table. The drawers can be ordered on the left (104L) or right (104R) of the patient at no additional cost. There are two hospital grade electrical outlets.

Choice of six upholstery and six panel colors (see back cover).

Specifications

Dimensions (top)	27" W, 54" L (68.6 x 137.2 cm)
Length with foot extension	72" (182.9 cm)
Height	32½" (82.6 cm)
Paper roll holder (hidden)	21" (53.3 cm)
Electrical outlets, duplex hospital grade	1 on left side
UL Listed	

Options

Pelvic Tilt (1°, 4°, 8°) and Drawer Heater
Stainless Steel Irrigation Pan (replaces treatment basin)
Urology Drain Pan (replaces treatment basin)

Accessories

Knee Crutches
Armboard
Extra Armboard Bracket
Procto Rest